



• Optional accessories shown

## Features

- FAST Chair (ordering number 60-4120-001) can be used as a stair chair, evacuation chair, or transport chair
- Ideal and cost effective solution utilizing proven Ferno EZ-Glide stair track system
- Stair track can be easily attached to or removed from the chair (minimizing weight to be carried by the operator)
- All-aluminium frame construction can carry patients of up to 250 kg weight
- Two 125 mm front swivel wheels, and two 150 mm rear fixed wheels with brakes
- Wheels have minimum debris collection while maintaining good traction and comfort, and provide very good maneuverability in confined spaces such as stair landings or tight corridors
- Head support and three restraints included
- Rear top handle is extendable and adjustable in four positions (1000 to 1580 mm)
- Front telescoping handles are adjustable in five positions (maximum extension of 340 mm)
- Soft vinyl cover is easy to clean or replace
- Vinyl cover can also be replaced with hard plastic seat and backrest (clip-on kit 65-4100-109)
- Optional Comfort Footbar, 4-point Restraint, Mesh Basket, Rear Support Handles, and Secure Lock Mount available
- FAST chair can be also ordered without stair-track (ordering number for a chair without stair-track is 60-4100-001), and the stair-track (manual or powered) can then be purchased separately later on if desired

## Specifications

Height min	Overall Width	Overall Length	Seat Width	Folded L,W, Thickness	Weight without track	Weight with track	Load Limit
1000 mm	516 mm	695 mm	415 mm	1000x516x215 mm	11,2 kg	15,0 kg	250 kg

### Accessories

FAST Chair Wall Mount - designed to securely hold the FAST Chair in both static or mobile application. Single hand release allows for quick chair fastening or detachment.

# 60-4190-001



4- Point Restraint for enhanced safety

# 62-4400-001



Hard Plastic Seat and Backrest

# 65-4100-109 (clip-on kit)



Rear Support Handles provide additional gripping points for the operator when lifting the chair or maneuvering over obstacles.

# 65-4122-001



Foldable Comfort Footbar increases patient's comfort and sense of security.

# 65-4130-001




Head Support – included free of charge

# 65-4126-001



Foot Restraint– included free of charge

# 65-0209-002

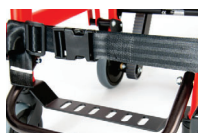

Foldable armrests for enhanced comfort.

# 65-4128-001



KIT (2PCS) ANGLED REAR HANDLES (adding 70mm extra space in-between) with folding-lock function # 65-4157-001


KIT (2PCS) ANGLED REAR HANDLES (adding 70mm extra space in-between) without folding-lock function # 65-4159-001

KIT (2PCS) REPLACEMENT ANGLED REAR HANDLES # 65-4158-001



Mesh Basket for storing accessories such as oxygen bottle.

# 65-4123-001

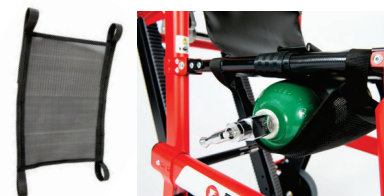

Manual stair-track for a chair purchased without it.

# 65-0402-001



FAST Big Wheels - suitable for off-road use, can be quickly and easily attached or removed and stored

# 65-0407-001-P



Powertraxx FAST - intelligent battery powered stair-track system that glides up to 227kg loads up or down the stairs without carrying or lifting, with up to 20 flights of stairs per charge. Chair is easy to operate when climbing or descending stairs, and you are always in control with intuitive and easy-to-use paddle-style controls.

# 65-0403-001



FAST chair has been designed to comply with and exceed international standards such as EN 1865-4:2012.

