



• *Optional accessories shown*

Features

- FST Chair (ordering number 60-4110-001) can be used as a stair chair, evacuation chair or transport chair
- Vinyl cover can also be replaced with hard plastic seat and backrest (clip-on kit 65-4100-110)
- Ideal and cost effective solution with Ferno EZ Glide track system
- All-aluminium frame construction can carry patients of up to 250 kg weight, and integrated flip-out stair-track allows gliding down the stairs without carrying or lifting
- Two 125 mm front swivel wheels, and two 150 mm rear fixed wheels with brakes
- Wheels have minimum debris collection while maintaining good traction and comfort, and provide very good maneuverability in confined spaces such as stair landings or tight corridors
- Head support, three restraints and comfort bar included
- Rear top handle is extendable and adjustable in four positions (1000 to 1580 mm)
- Front telescoping handles adjustable in five positions (maximum extension of 340 mm)
- Vinyl cover is easy to clean or replace
- Optional Comfort Footbar, 4-point Restraint, Mesh Basket, Rear Support Handles, and Secure Lock Mount available

Specifications

Height min	Overall Width	Overall Length	Seat Width	Folded L,W, Thickness	Weight	Load Limit
1000 mm	516 mm	695 mm	415 mm	1000x516x200 mm	13,8 kg	250 kg

Accessories

FAST/FST Chair Wall Mount - designed to securely hold the FST Chair in both a static or a mobile application. Single hand release allows for quick chair fastening or detachment.

60-4190-001



4-Point Restraint for enhanced safety

62-4400-001



KIT PLASTIC SEAT FST CHAIR

65-4100-110



Rear Support Handles provide additional gripping points for the operator when lifting the chair or maneuvering over obstacles.

65-4122-001



Comfort Footbar increases patient's comfort and sense of security.

Included Free of Charge

65-4130-001



Head Support – Included Free of Charge

65-4126-001




Foldable armrests for enhanced comfort.

65-4128-001



Foot Restraint– Included Free of Charge

65-0209-002

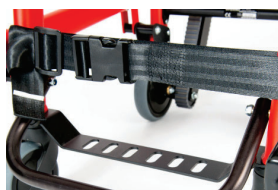

Mesh Basket for storing accessories such as oxygen bottle.

65-4123-001



KIT (2PCS) ANGLED REAR HANDLES (adding 70mm extra space in-between) with folding-lock function # 65-4157-001

KIT (2PCS) ANGLED REAR HANDLES (adding 70mm extra space in-between) without folding-lock function # 65-4159-001

KIT (2PCS) REPLACEMENT ANGLED REAR HANDLES # 65-4158-001



FST chair has been designed to comply with and exceed international standards such as EN 1865-4:2012.

