


• *Optional accessories shown*

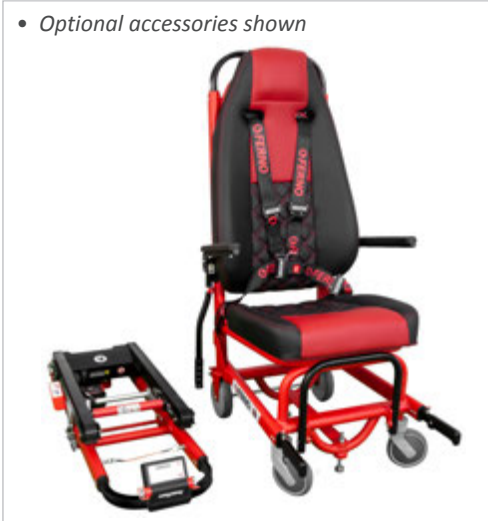

Features

- MAX Chair / non-foldable chair is used for transporting patients to the ambulance
 - Patients are transported in this chair in a seated position inside the ambulance
 - Frame constructed with welded aluminium for a long life and a carry capacity of up to 250 kg
 - Max Power is a ideal and cost effective solution of a removable Powertraxx with telescoping handle
 - That can be easily attached or removed from the chair (to minimize the weight that needs to be carried by the operator)
 - Red powder-coating finish prevents oxidation for clean hands and clothing
 - Anatomically-shaped seat
 - Easy to clean 100% Polyester, Vinyl cover can be disinfected
 - Armrest kit with anti-slip plastic coating
 - Two 125 mm front swivel wheels and two 150 mm rear fixed wheels with locks
 - Telescoping carrying handles at front and rear ends
 - Front telescoping handles adjustable in five positions (maximum extension of 340 mm)
 - Folding footrest
 - Optional MAX Power, IV Pole
- 60-4500-001 MAX CHAIR
 - 65-0409-001 MAX POWER
 - 60-0251-008 MAX FASTENER

Specifications

Overall Height	Overall Width	Overall Length	Seat Width	Weight without track	Load Limit
1225 mm	595 mm	830 mm	445 mm	20 kg	250 kg

Accessories

MAX POWER includes -

a telescoping handle and an intelligent stair-track battery powered system that easily glides up and down stairs with loads up to 227 kg without the need to carry or lift (up to 20 flights of stairs per charge). The chair is easy to operate when climbing or descending stairs with intuitive, easy-to-use paddle-style controls to keep you in control

Weight: 18kg

65-0409-001



IV Pole

This rugged and lightweight IV Pole offers an integrated height adjustment which allows operators to set their desired height for treatment

65-1500-031

Fastener

MAX Fastener

- tested and compliant with EN 1789 Standards
- stainless steel housing components
- automatic locking when loading chair into fastener
- manual release button (with mechanical release)

60-0251-008



MAX chair has been designed to comply with and exceed international standards such as EN 1865-1:2010+A1:2015 and EN 1789:2020.

